



Barwon Lift Recliner



Barwon Lift Recliner

The Barwon is a Lift-Recline Chair, the first choice for value in a high-quality Australian made lift chair. Choose from four standard chair size and an extensive range of feature adjustments to achieve a completely personalised fit and ensure comfortable usage. Every chair is fitted with carefully tested and selected electrical components and can be upholstered in your choice of premium fabric from Oscar's very own Versotela range.



FEATURES

- 130kg weight capacity
- Single motor wall saver
- Controller right side
- 9 volt battery
- 50mm casters
- Comes standard with Velcro attached back cushions
- Emergency battery backup
- Made to measure
- Available in wide range of Versotela fabrics
- 2 Year Warranty on mechanism, electrics and fabric
- 5 Year Warranty on timber frame
- 100% manufactured in Australia by Oscar Furniture




DESIGN

The Barwon comes standard with Velcro attached back cushions. The Velcro makes it easy to clean or replace cushions on site, and can help with minor adjustments to the cushion's position and fill.

Barwon three-roll back.


The top two rolls are filled with a fibrous filling that has a wonderful soft feel. Whilst the bottom roll is filled with a supportive foam. Barwon is the ultimate in softness and luxury.

ACTIONS

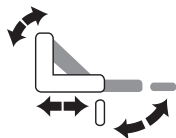

Single Motor (Wall-saver, Lift Assist)

The single motor wall saver operation of the Barwon allows the chair to move forward while the backrest reclines. This Oscar feature reduces the space that the chair occupies in the room because the reclined backrest does not hit the wall. The Lift, recline and footrest operate from one motor.

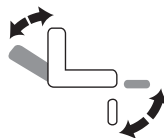

Dual Motor (Independent Recline)

Dual motor option lets the user recline the back without effecting the footrest, allowing a greater range of comfort.

SPECIFICATIONS



Seat Height

Seat Width

Seat Depth

Arm Height

Backrest Height

Scapular Height

Footrest Length

MINI SIZE

MIN ADJUSTMENTS	-30mm	-150mm	-150mm	-60mm	-180mm	Linked to back height	-50mm
MEASUREMENT	456mm	540mm	455mm	199mm	636mm	429mm	N/A
MAX ADJUSTMENTS	+20mm	+50mm	+75mm	+60mm	+150mm	Linked to back height	+50mm

SIZE A

MIN ADJUSTMENTS	-30mm	-150mm	-150mm	-60mm	-180mm	Linked to back height	-50mm
MEASUREMENT	476mm	540mm	490mm	180mm	647mm	440mm	N/A
MAX ADJUSTMENTS	+20mm	+50mm	+75mm	+60mm	+150mm	Linked to back height	+50mm

SIZE B

MIN ADJUSTMENTS	-30mm	-150mm	-150mm	-60mm	-180mm	Linked to back height	-50mm
MEASUREMENT	496mm	540mm	520mm	190mm	688mm	481mm	N/A
MAX ADJUSTMENTS	+25mm	+50mm	+50mm	+60mm	+150mm	Linked to back height	+50mm

SIZE C

MIN ADJUSTMENTS	-30mm	-150mm	-150mm	-60mm	-180mm	Linked to back height	-100mm
MEASUREMENT	546mm	540mm	545mm	190mm	718mm	511mm	N/A
MAX ADJUSTMENTS	+50mm	+50mm	+25mm	+60mm	+120mm	Linked to back height	+50mm

MODIFICATIONS

MEMORY FOAM

180KG WEIGHT CAPACITY

DUAL ACTION

SEAT WIDTH

SEAT HEIGHT

BACK HEIGHT

SEAT DEPTH

HEAD AND ARM REST COVERS

CONTROLLER LEFT SIDE

GLUTERAL RECESS

FOOTREST LENGTH

24 VOLT BATTERY

CUSHION CUTOUT

NO CASTERS

NO SHOWWOOD

INCREASE SIZE

TABLET ARM

ARM HEIGHT

MANUAL MECHANISM

LEG REST SIDE WALLS

GLUTERAL RECESS

REMOVABLE HEADREST



FABRIC SELECTION




VersoTela
Fabrics

All VersoTela fabrics come standard with a suite of high-end specifications, such as, antimicrobial coating, mold resistance, scratch resistance and fire resistance.

Our wide range of collections each boast unique properties that are carefully designed to suit every customer's needs. Whether it be built for use in a care facility, or as a functional but stylish addition to the living room, our VersoTela fabric options are sure to tick all the boxes when it comes to crafting your perfect chair.

Please get in touch with one of our friendly sales representatives for fabric samples and to find out more about which fabric is right for you.

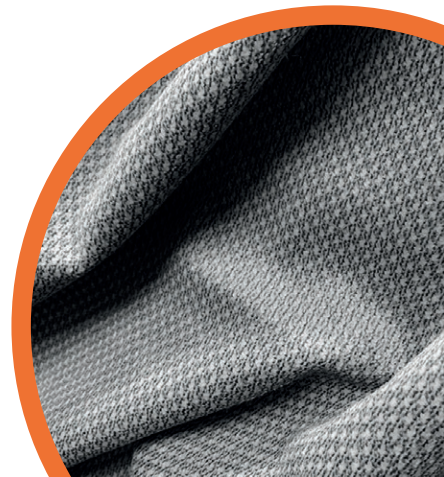

Optional Extras

Enhance Your Lift Chair Journey

COMFORT PACK

Elevate your Barwon chair to the pinnacle of comfort with our exclusive Comfort Pack upgrade. Experience a new level of personalised luxury, complete with integrated enhancements. Crafted with precision in our facility in Horsham, Australia, this add-on module is encased in a specially designed rigid foam unit for durability and reliability.



KEY FEATURES:

- **Forward Headrest Tilt:** The headrest actuator's low pivot point is ingeniously designed to mimic natural head and neck movements, ensuring optimal support and comfort.
- **Lumbar Adjustment:** Enjoy enhanced lumbar support thanks to the integrated air bladder and pump unit. Activate the pump, and it silently inflates to provide the exact level of support you desire. This technology, initially pioneered in the automotive industry for car seat lumbar adjustment, remains at the forefront of innovation in this field.
- **Convenient Built-in Controller:** Effortlessly operate these advanced functions with the conveniently located controller on the right-hand side of the chair. (Lift and recline functions remain controlled by traditional hand control).
- **Service-Friendly Design:** The Comfort Pack is designed with easy servicing in mind, aligning with Oscar Group's eco-friendly focus on durable and repairable products.



CHAIR COVERS

Discover the perfect blend of style and functionality with our Barwon chair covers, meticulously crafted to offer a bespoke touch while safeguarding your cherished lift chair from spills and pet mishaps.



KEY FEATURES:

- **Universally Designed:** Our chair covers are designed to flawlessly fit Barwon Mini, A, B, and C chairs, ensuring a snug and tailored appearance.
- **Quality Fabric Options:** Choose from our extensive range of Versotela fabrics, allowing you to match your chair cover with your interior decor effortlessly.
- **Customisable Design:** Tailor your chair's protection to your exact needs. Select between a full cover, as elegantly displayed, or opt for a back seat and leg cover without arms.



Contact us today:
sales@oscarfurniture.com.au
03 5381 1404

6 King Drive,
Horsham VIC 3400

www.oscarfurniture.com.au

